



Stealth Products[®]



All-Round Lite Joystick

Owner's Manual
mo-Vis All-Round Lite Joystick for R-Net


moving forward together

Table of Contents

| | | | |
|--|------------|--|-----------|
| 1.0 Customer Satisfaction | i | 9.0 Installation Instructions | 7 |
| 1.1 General..... | i | 9.1 Preparations..... | 7 |
| 2.0 Important Information..... | ii | 9.2 Tools | 7 |
| 2.1 Important Information | ii | 9.3 Installation Plan | 8 |
| 2.2 Disposal | ii | 9.4 Joystick Installation | 8 |
| 2.3 Supplier Reference | ii | 9.5 Satellite Twister Install..... | 9 |
| 3.0 Introduction | iii | 9.6 Connection Cable Install | 10 |
| 3.1 Ordering Documentation | iii | 10.0 Configuration Software | 11 |
| 4.0 Warranty | iv | 10.1 Software Capabilities | 11 |
| 4.1 In Case of Product Failure..... | iv | 10.2 Software Download..... | 11 |
| 5.0 Warning Labels..... | v | 11.0 R-Net Setup | 12 |
| 5.1 Warning Labels | v | 11.1 R-net Setup..... | 12 |
| 5.2 Limited Liability..... | v | 11.2 R-net Further Setup..... | 12 |
| 5.3 Testing | v | 11.3 R-net Para. Settings..... | 14 |
| 6.0 Design and Function..... | 1 | 12.0 Software Configuration | 17 |
| 6.1 Intended Use | 1 | 12.1 mo-Vis Config. Para. Settings..... | 17 |
| 6.2 Features | 1 | 13.0 Testing..... | 19 |
| 6.3 Mounting | 1 | 13.1 Testing..... | 19 |
| 7.0 Parts and Accessories | 2 | 13.0 Error Codes and Fault Logs | 21 |
| 7.1 Micro Joystick Package..... | 2 | 13.1 Error Codes | 21 |
| 7.2 Micro Joystick Accessories | 2 | 13.2 Fault Logs | 22 |
| 7.3 Micro Joystick Prot. Covers | 3 | 15.0 Technical Data | 24 |
| 8.0 Functioning..... | 4 | 15.1 Technical Data..... | 24 |
| 8.1 Joystick Operation..... | 4 | 16.0 First Time Use | 25 |
| 8.2 Joystick Safety | 5 | 16.1 Dealer Assistance..... | 25 |
| 8.3 Connection Options..... | 6 | 16.2 User Testing | 25 |
| 8.4 EM Compatibility..... | 6 | 16.3 Conditions of Use..... | 25 |
| | | 17.0 Maintenance..... | 26 |
| | | 17.1 Maintenance..... | 26 |
| | | 17.2 Cleaning | 26 |
| | | 17.3 Monthly Check | 26 |
| | | 18.0 Notes | 27 |

Stealth Products is committed to 100% customer satisfaction. Your complete satisfaction is important. Please contact us with feedback or to suggest changes that may help us improve the quality and usability of our products.

You may reach us at:



Stealth Products

104 John Kelly Drive, Burnet, TX 78611


Phone: (512) 715-9995

Toll Free (800) 965-9229

Fax: (512) 715-9954

Toll Free (800) 806-1225

info@stealthproducts.com

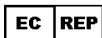

MDSS GmbH

Schiffgraben 41

30175 Hannover, Germany



General

1.1

Read and understand all instructions prior to the use of this product. Failure to adhere to instructions and warnings in this document may result in property damage, injury, or death. Product misuse due to failure to adhere to the following instructions will void the warranty.

Immediately discontinue use if any function is compromised, if parts are missing or loose, or if any component shows signs of excessive wear. Consult your supplier for repairs, adjustments, or replacements.

Important Information

2.1

All persons responsible for fitting, adjustment, and daily use of the devices discussed in these instructions must be familiar with and understand all safety aspects of the devices mentioned. In order for our products to be used successfully, you must:

- read and understand all instructions and warnings;
- maintain our products according to our instructions on care and maintenance; and,
- ensure devices are installed and adjusted by a trained technician.

Disposal

2.2

For scrapping, please adhere to your local waste legislation. Dispose of obsolete electronic parts responsibly and in accordance with local recycling regulations.

Supplier Reference

2.3

Supplier:

Telephone:

Address:

Purchase Date:

Model:

The installation instructions will guide you through the options and possibilities with the product. The **All-Round Lite Joystick** should be adjusted with the mo-Vis Configurator Software.

Instructions are written with the expressed intent of use with standard configurations. They also contain important safety and maintenance information, as well as describe possible problems that can arise during use. For further assistance or more advanced applications, contact your supplier or **Stealth Products** at (512) 715-9995 or toll free at (800) 965-9229.

Always keep the operating instructions in a safe place so they may be referenced as necessary.

All information, pictures, illustrations, and specifications are based on the product information that was available at the time of printing. Pictures and illustrations shown in these instructions are representative examples and are not intended to be exact depictions of the various parts of the product.

**CAUTION**

This product is designed to be fitted, applied, and installed exclusively by a healthcare professional trained for these purposes. The fitting, application, and/or installation by a non-qualified individual could result in serious injury.

Ordering Documentation

3.1

Additional copies of this manual may be downloaded by accessing the **Stealth Products** website:

<https://stlpro.site/stealth-docs>

and searching

"mo-Vis All-Round Joystick"

in the search bar at the top of the page.

Our products are designed, manufactured, and produced to the highest of standards. If any defect in material or workmanship is found, **Stealth Products, LLC** will repair or replace the product at our discretion. Any implied warranty, including the implied warranties of merchantability and fitness for a particular purpose, shall not extend beyond the duration of this warranty. **Stealth Products, LLC** does not warrant damage due to, but not limited to:

- misuse, abuse, or misapplication of products; and/or
- modification of product without written approval from **Stealth Products, LLC**.

Any lack or alteration of serial number, where applicable, will automatically void this warranty.

Stealth Products, LLC is liable for replacement parts only. **Stealth Products, LLC** is not liable for any incurred labor costs.

No person is authorized to alter, extend, or waive the warranties of **Stealth Products, LLC**.

Covers/Soft Goods: 2 years

Hardware: 5 years

Electronics: 3 years

In Case of Product Failure

4.1

In the event of product failure covered by our warranty, please follow the procedure outlined below:




1. Call **Stealth Products** at (512) 715-9995 or toll free at (800) 965-9229.
2. Request a Return Authorization form (RA) from the Returns Department and follow documentation instructions.

IN THE CASE OF ANY SERIOUS INCIDENT OCCURRING IN RELATION TO THIS PRODUCT, A REPORT TO STEALTH PRODUCTS, LLC SHOULD BE MADE IMMEDIATELY.

Warning Labels

5.1

Warnings are included for the safety of the user, client, operator, and property. Please read and understand what the signal words **DANGER**, **WARNING**, **CAUTION**, **NOTICE**, and **SAFETY** mean, and how they could affect the user, those around the user, and property.

| | |
|---|--|
|  DANGER | Identifies an imminent situation which, if not avoided, may result in severe injury, death, and property damage . |
|  WARNING | Identifies a potential situation which, if not avoided, may result in severe injury, death, and property damage . |
|  CAUTION | Identifies a potential situation which, if not avoided, may result in minor to moderate injury and property damage . |
| NOTICE | Identifies important information not related to injury, but possible property damage . |
| SAFETY | Indicates steps or instructions for safe practices, reminders of safe procedures, or important safety equipment that may be necessary. |

Limited Liability

5.2

Stealth Products, LLC accepts no liability for personal injury or damage to property that may arise from the failure of the user or other persons to follow the recommendations, warnings, and instructions in this manual.

Stealth Products LLC does not hold responsibility for final integration of final assembly of product to end user.

Stealth Products, LLC is not liable for user death or injury.

Testing

5.3

Initial setup and driving should be done in an open area free of obstacles until the user is fully capable of driving safely.

Intended Use

6.1

The **All-Round Lite Joystick** is a compact joystick with half the normal force and one-third the normal throw of the standard All-Round joystick. Both joysticks are compact joysticks that require more force than the Micro or Multi joysticks, and both are used as proportional joysticks able to be controlled by hand, chin, or by an attendant.

An R-net Joystick has a communication cable with an R-net connector. The cable can be directly plugged into the electronics of the wheelchair. The joystick then controls the wheelchair in all its functions (driving, positioning, menu navigation, etc.).

Features

6.2

mo-Vis' **All-Round Lite Joystick** features:

- full proportionality and a compact design;
- a light amount of throw and required force
- the ability to be used as an attendant joystick, or operated by hand or chin;
- the ability to link with either one or two switches;
- full adjustability; and,
- easy configuration via the free mo-Vis Configuration Software.

Mounting

6.3

The **All-Round Lite Joystick** can be mounted at any location on the wheelchair. By attaching the Power Chin Boom to the chair, the All-Round Lite can be used as a chin joystick.*

The **All-Round Lite Joystick** is compatible with Stealth joystick mounts, including:

- the swing-away SM-660M;
- The GAT400 series GAT472P-M; and,
- The ARMS260 series mounts.

*** Note:** A 6mm mounting set is required to mount the All-Round Lite to the Power Chin Boom.

All-Round Lite Joystick Package

7.1

The following is included in the **All-Round Lite Joystick** package:



| Product Description: | Part Number: |
|---|--------------|
| All-Round Lite Joystick Set, with integrated cabling and interface. | IDM-ARLITE-9 |

All-Round Lite Joystick Accessories

7.2

For mounting and further personalization, the following accessories are available (and are to be ordered separately):

| | Product Description: |
|---|---|
|  | All-Round/All-Round Lite joystick ball. |
|  | Satellite Twister (includes All-Round/All-Round Lite insert. (available in red, black, yellow, blue or green) |

All-Round Lite Joystick Protective Covers

7.3

The connections of the **All-Round Lite Joystick** unit have protective covers to avoid the accumulation of dust or moisture during use.

WITH PROTECTIVE COVERS IN PLACE:



WITH PROTECTIVE COVERS REMOVED:



When using the joystick, an LED on the All-Round Lite will indicate the unit is functioning properly. The LED will flash if an error occurs. (The LED also indicates the forward direction of the joystick.)

Please refer to the *Error Codes* section of this manual for a list of errors indicated by a blinking LED on the **All-Round Lite Joystick**.

NOTICE

Before inserting a connector into the joystick, remove the protective cover from the desired input port.



CAUTION

Keep the protective cover in its place when not using a specific port.

NOTICE

Not all items in the Multi Joystick will be used to mount the joystick onto a Stealth Product. All of these items are compatible with the mo-Vis Power Chin Boom.

Stealth Products can customize mounting hardware configurations in order to create the best possible fit for a client.

Joystick Operation

8.1

The movements of the joystick translate accordingly into the motion of the wheelchair or the navigation of its menu. The wheelchair is most commonly controlled by the joystick in the following manner:

- **Direction:** The joystick is pushed in the direction the wheelchair is intended to move, and the wheelchair moves in that direction.
- **Speed:** The further the joystick is pushed in any direction from its default (center) position, the faster the wheelchair moves in that direction.
- **Stop:** When the joystick is released from being pushed in a direction, it returns to its default (center) position and the wheelchair comes to a stop.



WARNING

When the LED flashes to indicate an error and/or after every incident with the wheelchair or the mo-Vis device, contact Stealth immediately to perform a function test.



WARNING

The **All-Round Lite Joystick** operates with high sensitivity, requiring sufficient movement coordination and control of force. When operating the joystick, take care to avoid erratic or unsteady movements.



WARNING

When operating the **All-Round Lite Joystick**, always follow the instructions provided by the wheelchair manufacturer regarding safely operating the chair on different surfaces and angles of ascent and descent.



WARNING

Users should not use the **All-Round Lite Joystick** as the only support for a hand or limb, as irregular movements may disrupt controls.



CAUTION

Avoid hitting obstacles during driving.

Joystick Safety

8.2

Do not use the **All-Round Lite Joystick** if any of the following conditions apply:

- The joystick handle is damaged, cracked, or missing.
- The joystick does not return to its default (center) position independently.
- The joystick does not move to and from its default (center) position smoothly.

In order to change the joystick tip, gently pull on the tip until it lifts off the joystick's metal rod. Place the replacement tip on the metal rod until it sits firmly in place.



WARNING

Always power off the device when the user is no longer in the wheelchair to avoid unauthorized access or accidental movements.

SAFE OPERATION

It is mandatory to have an on/off switch on the wheelchair which immediately shuts down the chair power and its electronics. It must be within easy reach of the user while operating the wheelchair. This allows the user to instantly stop the wheelchair in case of problems or an emergency.



WARNING

The on/off switch must be available to the user at all times while operating the wheelchair.



WARNING

Do not cover or block the device to avoid uncontrolled behavior of the device or wheelchair.

If the wheelchair responds in an unexpected way during use, the user should immediately release the joystick and/or use the power on/off switch.

Connection Options

8.3

The **All Round Lite Joystick** can connect one or two switches to either the joystick or the interface. One connection is to power the on/off function and the other is to function as an input or user switch for the wheelchair.

The interface has a cable with an R-net BUS connector compatible with the wheelchair electronics. The joystick then controls the wheelchair and all of its functions, including driving, positioning, menu navigation, etc. It is also capable of powering the wheelchair on or off, as well as issuing mode commands.

Jack on/off (power) = red
Jack mode (input) = yellow



Electromagnetic Compatibility

8.4

The electronics of a power wheelchair and its options can be affected by external electromagnetic fields (for example, from mobile telephones). Similarly, the electronics of the wheelchair or options themselves can also emit electromagnetic fields that can affect the immediate surroundings (for example, certain alarm systems in businesses).

This product complies with the limit values for Electromagnetic Compatibility (EMC) with respect to power wheelchairs, as set out in the harmonised standards for the EU in the Medical Device Regulation, No. 2017/745.

Preparations

9.1

Only a qualified service technician should install the **All Round Lite Joystick** and its accessories.



WARNING

An incorrect installation of wheelchair electronics may cause damage to the hardware and/or injury to the user.



CAUTION

Any connection must always be secured with all the delivered screws.
Only use the screws provided with the device.



CAUTION

Before you begin the installation: Please check the packaging to verify all items have been included. Make sure you have all necessary documentation and knowledge to complete the installation.
Check the condition of the device.

Tools

9.2

Use the proper tools to install and adjust the **All Round Lite Joystick** to the desired mounting hardware.



CAUTION

The use of improper tools may damage the device.
Do not over tighten the screws.
Doing so may damage the joystick unit.

Installation Plan

9.3

Set up an installation plan before beginning the installation. This plan should take into account:

- where the **All-Round Lite Joystick** will be placed;
- how the **All-Round Lite Joystick** will be operated (i.e. by hand, chin or attendant...); and,
- **All-Round Lite Joystick** parameter settings.



CAUTION

Any connection must always be secured with all the delivered screws.
Only use the screws provided with the device.

Joystick Installation

9.4

The **All-Round Lite Joystick** can be installed to previously ordered mounting hardware simply and easily. Installation is only a matter of connecting the cabling to the wheelchair electronics.

When installing the joystick, place the joystick in the desired location and secure all wires and cabling, ensuring that no cables or wires will obstruct the joystick during use. If necessary, before connecting the joystick to the wheelchair electronics, first attach a power on/off switch and/or a mode (in) switch.



CAUTION

Only install this product to wheelchairs onto which the chair manufacturer allows the installation of third-party accessories.



WARNING

Ensure all cabling is mounted in such a way to avoid rubbing against other equipment, causing unnecessary and excessive wear and tear.



CAUTION

Before inserting a connector, remove the input's protective cover.
Do not discard the protective cover.

If the **All-Round Lite Joystick** is not mounted in a position in which the USB connector is pointing toward the user, the joystick direction can be adjusted in increments of 90° with the mo-Vis Configurator Software.

- The software must be installed on a PC and ready to use.
- With Dealer-level access, several parameter settings can be changed, such as the mounting position. Dealer-level access requires a password, which can be easily obtained by contacting Stealth Products, LLC.

Satellite Twister Installation

9.5

Up to two Satellite Twisters can be mounted to an **All-Round Lite Joystick**. Satellite Twisters come with All-Round Lite inserts which must be used when mounting.

The steps to mounting a Twister to the joystick are as follows:

1. Determine the best location(s) for the Twister(s).
2. Unscrew and remove the cover for the Twister connection with a Phillips screwdriver.
3. Place the Twister firmly into the open slot.
4. Tighten the screw to secure the Twister to the joystick.
5. (If necessary, repeat these steps for the second Twister.)
6. Connect the cabling for the Twister to the joystick.



CAUTION

All wheelchair electronics must be switched off during installation.



CAUTION

The **All Round Lite** is a Class I Medical Device (MDR 2017/745). All accessories, including switches, must also be compliant with MDR 2017/745.



CAUTION


Excessive tightening of screws can cause damage to the joystick unit.

Connection Cable Installation


9.5

Insert the connector into its mate and push firmly. The connector is inserted correctly if the yellow pins are completely invisible. To disconnect the cable, firmly hold the connector housing and pull the connectors apart.


Disconnected

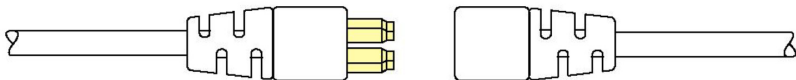

Correctly Connected

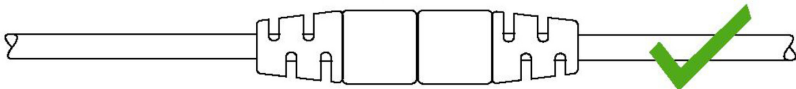

Incorrectly Connected

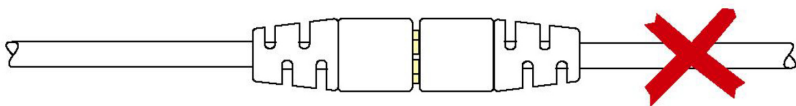

CAUTION

Before inserting a connector, remove the input's protective cover.
Do not discard the protective cover.

Software Capabilities

10.1

The **All-Round Lite Joystick** is adjustable with the mo-Vis Configurator Software.

- The mo-Vis Configurator Software is compatible with Windows 7, Windows 8, and Windows 10 operating systems.
- The mo-Vis Configurator Software is not compatible with Android or iOS systems.



WARNING

Changes in the parameter settings may cause damage to the device or power chair, or may cause injury to people.



CAUTION

Always change parameters and test the outcome without anyone sitting in the chair,.

Software Download

10.2

The mo-Vis Configurator Software can be downloaded on Stealth Products' website (www.stealthproducts.com). For details on how to install and use the software, the user is advised to consult the software manual, included alongside the download.

Please contact your dealer for information about downloading the R-net PC Programming software.

R-net Setup

11.1

For mo-Vis R-net joysticks to operate properly, some R-net parameters need to be set. For information about R-net parameters, please refer to the Curtiss-Wright manual *SK77981-14 | R-net Technical Manual | Chapter 3 -- Programming*.

- *R-net Profile Management > Input Device Type (Ch. 3 - 4.4)*: In most cases, this parameter needs to be set to *Universal* for all mo-Vis joysticks. In combination with certain other devices, it is advised to set this parameter to *JSM*. Please refer to the installation manual of the other devices for more information.

- *R-net Controls > Global > Profile Button (Ch. 3 - 7.5)*: The default setting for this parameter is *Profiles*. To access the control modes for other wheelchair functions, such as seating positions, select *Profiles/ Modes*.

R-net Further Setup

11.2

With the R-net PC programming, you can fine-tune joystick performance and functionalities such as:

- Joystick Throw

JOYSTICK THROW

Throw, as it relates to joystick operation, is the position of a joystick, relative to its default (center) position, when it is stopped in a particular direction. *Throw distance*, then, is the maximum distance a joystick shaft can be moved in any particular direction.

As it relates to programming, throw will correlate to the speed at which the wheelchair will move in a given direction; a joystick thrown maximally in a direction will cause the chair to move in that direction at its maximum speed.

When calibrating a joystick, its throw should be set to the maximum distance the user is capable of moving the joystick shaft so that the user's full efforts correspond with the chair's maximum speed output, i.e. a setting of 50% will correspond to full throttle at half the joystick's throw distance.

CHANGING THE THROW SETTINGS

1. Plug the R-net programming dongle into an available port in the R-net system.
2. Power on the chair.
3. Connect the PC to the Programming Dongle with the USB cable.
4. Start the R-net PC Programming Software on your PC.
5. Select *Read from Controller* to load the information from the R-net system to the Configurator software.
6. To change the setting for forward throw, select *Controls > Joystick > Joystick Forward Throw*.
7. Enter the value for the parameter. (Value range is 25% to 100%, in increments of 1%)
8. Repeat steps 6 and 7 to edit the settings for *Joystick Reverse Throw*, *Joystick Left Throw*, and *Joystick Right Throw*, as necessary.
9. Select *Write to Controller* to save the changes to the R-net system.
10. Power cycle the chair (turn it off, then back on) to load the new parameters.

R-net Parameter Settings

11.3

With the mo-Vis Configurator software you can change the parameters of the Micro Joystick. In the dealer level, you will be able to change a number of parameters.

| Global parameter | Supported | Version | Solutions |
|---------------------------------|-----------|---------|-----------|
| Momentary screens enabled | N | | |
| Change profile while driving | N | | |
| Change speed while driving | N | | |
| Actuator switches while driving | N | | |
| Speed adjust | N | | |
| Profile button | Y | V02.00 | |
| Actuator endstop beep | N | | |
| Sounder volume | Y | V02.03 | |
| Start-up beep | N | | |
| Lock function enabled | N | | |
| Reverse driving alarm | Y | V02.00 | |

| Global parameter | Supported | Version | Solutions |
|--------------------------------|-----------|---------|-----------|
| Emergency stop switch | N | | |
| OBP keycode entry | N | | |
| Power-up mode | N | | |
| External profile jack function | N | | |
| Profile/mode Jack detect | N | | |
| On/off jack detect | N | | |

| Profile parameter | Supported | Version | Solutions |
|------------------------|-----------|---------|---|
| Joystick forward throw | Y | V02.01 | |
| Joystick reverse throw | Y | V02.01 | |
| Joystick left throw | Y | V02.01 | |
| Joystick right throw | Y | V02.01 | |
| Joystick deadband | N | | Use deadband setting in the mo-Vis Configurator |

| Profile parameter | Supported | Version | Solutions |
|---------------------------|-----------|---------|--|
| Invert left/right JS axis | Y | V02.01 | |
| Invert fwd/rev JS axis | Y | V02.01 | |
| Swap joystick axis | N | | Use rotate function in the mo-Vis Configurator |
| Change mode while driving | Y | V02.00 | |
| Sleep timer | Y | V02.00 | |
| Standby timer | Y | V02.00 | |
| Switch to standby | Y | V02.00 | |
| Mode selection in standby | Y | V02.01 | |
| Standby in modes | Y | V02.01 | |
| Standby forward | Y | V02.00 | |
| Standby reverse | Y | V02.00 | |
| Standby left | Y | V02.00 | |
| Standby right | Y | V02.00 | |
| Remote selection | Y | V02.03 | |
| Background | N | | No screen to support this. |

mo-Vis Configurator Parameter Settings

12.1

The mo-Vis Configurator Software allows users to change the parameters of the **All Round Lite Joystick**, depending on the level of profile assigned.

MOUNTING SETTINGS (DEALER LEVEL)

| Setting | Description | Parameters | |
|---------------------------|--|-------------|------|
| <i>Mounting Direction</i> | Changes the joystick axis in 90-degree increments clockwise, providing compensation for joystick rotation. | Default: | 0° |
| | | Min.: | 0° |
| | | Max.: | 270° |
| | | Increments: | 90° |

DEADBAND SETTINGS

| Setting | Description | Parameters | |
|-----------------|--|------------|-----|
| <i>Deadband</i> | Neutral zone around the joystick's default (center) position | Default: | 10% |
| | | Min.: | 10% |
| | | Max.: | 50% |

(Note: Most wheelchair electronics have their own default deadband setting. It is recommended to keep the deadband setting at its default in the mo-Vis Configurator.)

TILT

The **All Round Lite Joystick** unit is able to detect tilt (in reference to the gravitation of Earth). The wheelchair will stop driving when the *Enable* parameter is set to *On* and the joystick is tilted more than 70° in any direction. The joystick becomes active again if the angle of tilt falls below 60°. The use of this functionality makes sense when a mo-Vis joystick is integrated into a tray or other tilting surface or mounting aid.

(Note: the default setting for the tilt sensor is *Off*. The tilt sensor response time is 1s. or less.)

**CAUTION**

An incorrect programming of the wheelchair or device electronics may cause damage to the devices or injury to the user.

COMPENSATION MODE

| Setting | Description | Parameters | |
|--------------------------|--|---------------|------------------------------------|
| <i>Compensation Mode</i> | Defines the activation of the compensation algorithm | <i>Off</i> | Always deactivated |
| | | <i>Manual</i> | (De)activated manually by the user |
| | | <i>On</i> | Always activated |

(Note: *Front wheel-driven wheelchairs are less affected by this factor than rear wheel-driven chairs.*)

| Setting | Description | Parameters | |
|-----------------------|--|----------------|-------------------|
| <i>X Compensation</i> | Limits the steering reaction if the wheelchair overreacts on steering commands. (This may happen due to compensation in the X direction that prevents the wheelchair from overcoming an obstacle due to torque loss.) | <i>None</i> | No Compensation |
| | | <i>1/2 Y</i> | X=50% of Y level |
| | | <i>2/3 Y</i> | X=66% of Y level |
| | | <i>Equal Y</i> | X=100% of Y level |

(Note: *The values are based on the compensation in the Y direction.*)

| Setting | Description | Parameters | |
|-------------------------------|---|---------------|------------------------|
| <i>(De)activation Pattern</i> | Active if Compensation Mode is set to Manual (to activate/deactivate, nudge the joystick forward, then backward, then forward. | <i>Fast</i> | Fast execution speed |
| | | <i>Normal</i> | Normal execution speed |

(Note: *Prevents the wheelchair from creeping after the joystick is released.*)

COMPENSATION FACTOR

| Setting | Description | Parameters | |
|----------------------------|---|------------------|-------------------------------------|
| <i>Compensation Factor</i> | Lowers forward/backward (Y) driving speed on rough terrain. (Wheelchair type and weight impact the compensation factor.) | <i>Very Weak</i> | Almost no deceleration |
| | | <i>Weak</i> | Decelerates less than <i>Normal</i> |
| | | <i>Normal</i> | Default Setting |
| | | <i>Strong</i> | Decelerates more than <i>Normal</i> |

(Note: *If the joystick is mounted in a location susceptible to bumps and shocks, such as a wheelchair arm, the compensation setting should be set to 'Weak' or 'Very Weak'.*



WARNING

After adjusting the joystick settings, always test the operation and functioning of the wheelchair before use.

Testing

13.1

After installing the **All-Round Lite Joystick**, execute the following tests, in order, before the wheelchair is delivered or put into service.

- a visual check to ensure that the joystick is intact;
- an operational test;
- a driving test; and,
- a test of the chair's stopping ability.

VISUAL TEST

Check to ensure:

- the joystick is not bent or damaged;
- none of the joystick's housing, cabling, or connectors are damaged; and,
- the joystick returns to its default (center) position upon release.



WARNING

The device should always be tested without any person sitting in the wheelchair after every alteration of the physical installation or adjustment of the parameters.

OPERATIONAL TEST

To perform an operational test, complete the following steps:

1. Activate the wheelchair operating system.
2. Check for any error messages.
3. Gradually push the joystick forward until you hear the chair's parking brakes release (The wheelchair may start moving.).
4. Immediately release the joystick. You should hear the parking brake activate within a few seconds.
5. Repeat steps 3 and 4 three times, each time moving the joystick forward, left, right, and backward.
6. Check to ensure the power switch and the mode switch are functioning properly, and within easy reach of the end user.

DRIVING TEST

Test drive the wheelchair, ensuring that:

- the chair and all its operations are fully functional in all positions the user may use the **All-Round Lite Joystick** and any connected switches;
- no cabling or connectors are positioned in such a way that they may become damaged during operation; and,
- the joystick returns to its default (center) position upon release.



WARNING

Only perform an operational test on a level surface, in an open area free of obstacles.

STOPPING TEST

Drive the wheelchair at forward at full speed, then use the power on/off switch to shut the wheelchair down.

- The wheelchair should not stop abruptly, but should gradually slow to a stop.

Error Codes

14.1

The illumination of the LED on the joystick interface indicates the operational status of the joystick. The table below describes the states the LED can indicate:

| R-net Status | Tilt Parameter | Joystick Tilted | LED Status |
|--------------|----------------|-----------------|--------------------------|
| Out of Focus | - | - | LED Flashes as heartbeat |
| In Focus | Enabled | N | LED on |
| In Focus | Enabled | Y | LED Flashes as heartbeat |
| In Focus | Disabled | Y | LED on |
| In Focus | Disabled | N | LED On |
| Configuring | - | - | LED flashes quickly |
| Power Cycle | - | - | LED flashes as heartbeat |
| Error | - | - | LED flashes error code |

When a fault occurs, the LED of the **All Round Lite Joystick** will start to flash. A long delay is followed by a number of flashes with a short delay. Count the number of flashes and look up the corresponding error message in the table below:

| Flash Count: | Reason: | Required Action: |
|--------------|--|--|
| 1 | - | - |
| 2 | Connection Cable/Driver Issues | Check/Replace Cable to Wheelchair and/or check/replace sensor cable (if available), or Replace PCB |
| 3 | Power Supply | Check R-net Cable and/or Replace PCB |
| 4 | Joystick/Sensor Fault | Check/Replace Joystick |
| 5 | - | - |
| 6 | ADC - Internal Analog to Digital Converter | Replace PCB |
| 7 | Test Flag Failed or Diagnostic Failed | Redo Tests and/or Replace PCB |
| 8 | CPU Fault | Replace PCB |
| 9 | Scheduler Fault | Update Software or Replace PCB |
| 10 | Coding Error | Update Software or Replace PCB |

Fault Logs

14.2

FAULTS

The configurator software maintains a fault log with fault counters. Each time a specific fault occurs, its counter will increase by one. (The maximum count value is 254.) The fault log can be accessed by the configurator (which requires Dealer-level profile access), which is then able to clear one or all of the fault counters. Reference the information presented in the following chart:

| Fault: | Reason for Fault: | Required Action: |
|-------------------------|---|--|
| CPU Error: RAM | CPU consistency check failed | Replace PCB. |
| CPU Error: FLASH | | |
| CPU Error: EEPROM | | |
| Run Error: Scheduler | Firmware consistency check failed | <ul style="list-style-type: none"> - Update software. - Replace PCB. |
| Code Error: Framework | | |
| Code Error: Application | | |
| MSP Command Corrupt | Corrupt Command Received | Connection with PC (Configuration Software) went wrong, try again. |
| MSP Command Unknown | Unknown Command Received | Connection with PC (Configuration Software) went wrong. Update Firmware or Update Configurator Software. Try again. |
| MSP Sub Command Unknown | Unknown Sub Command Received | |
| MSP Argument Invalid | Invalid Argument Received | |
| MSP Device Not Ready | Device not ready to receive MSP command | |
| MSP Device Wrong State | Device unable to receive MSP command in current state | |
| PCB Test Failed | Factory test failed | A fault occurred during factory testing. |
| Assembly Test Failed | | |

| Fault: | Reason for Fault: | Required Action: |
|-----------------------------|---|--|
| Field Test Failed | Field test failed (Calibration). | A fault occurred during Field Testing (calibration). |
| Test Flag Check | One or More Test flags Not Set. | Redo test(s) and/or Replace PCB. |
| ADC | ADC Conversion error | - Check R-net Cable, Replace PCB Interface. |
| R-net UART Overflow | UART send queue is full. | Replace PCB. |
| R-net UART Underflow | UART receive queue is empty. | |
| R-net Communication Timeout | Maximum number of packet re-transmissions is reached. | |
| R-net Tx Overflow | Packet Transmit Buffer is Full. | |
| R-net Rx Overflow | Packet Receive Buffer is Full | |
| R-net Invalid Seq Nr | Received a Packet with an Invalid Sequence Number | |
| R-net Data Packet Error | Data Packet ACK mismatch | |
| R-net Descr Error | Invalid Packet Data Descriptor | |
| R-net API Version Error | The R-net Chip Contains an Invalid API Version | Load the latest API version into the R-net chipset with the R-net Dongle. |
| R-net Chip Tripped Error | The R-net Chip has encountered an Internal Error. | Internal chipset error: Replace PCB R-net system error: solve R-net system error. |
| Communication | Communication with the sensor (joystick failed) | - Check Sensor (joystick) cable. - Replace sensor (joystick). |
| Joystick (Sensor) | Faulty sensor (joystick) | |
| Accelerometer | Accelerometer failed | Replace PCB Interface. |

Technical Data

15.1

JOYSTICK CONNECTORS

- 1/8" (3.5mm) mono jack input (on/off)
- 1/8" (3.5mm) mono jack input (mode)
- Mini USB
- Rnet BUS cable

JOYSTICK DIMENSIONS

- 96.5 x 66 x 60 mm / 3.8 x 2.60 x 2.36 in. (H x W x D)

JOYSTICK DATA

- Operational Force: 120gf (.26 lbf)
- Angle from Center to End: 8°
- Length of Joystick from Mechanical Pivot: 1.96 in. (50mm)
- Travel from Center: .27 in. (7mm)
- Maximum Vertical Load: 30591g (67.44 lb.)
- Maximum Horizontal Load: 15295g (33.72 lb)
- Expected Life: over 2,800,000 rotations
- Mass: 315g
- Interface Cord Length: 70.9 in. (180cm)
- Operating Temperature: -13°F to 122°F (-25°C to 50°C)
- Storage Temperature: -40° to 149°F (-40°C to 65°C)
- Immunity Level (ISO7176-21): 20V/m (26MHz to 2.5GHz)
- Emission Level (ISO7176-21): CISPR 11, Class B
- ESD (ISO7176-21): 8kV Air (4kV Contact)

Dealer Assistance

16.1

During first time use by the client, it is advised that the dealer or service technician not only install the device, but also explains the configuration and different possibilities to the customer (i.e., the user and/or the attendant). If needed, the dealer can make final adjustments.

User Testing

16.2

It is important that the customer is fully aware of the installation of the **All-Round Lite Joystick**, how to use it, and how it can be adjusted to fit the client comfortably. As a dealer, proceed as follows:

- Explain and show the customer how you have executed the installation and explain the function.

Have the user test the position of the **All-Round Lite Joystick**.

- Is the hardware in the proper position for the client?
- Can the user safely operate all controls with minimal effort?

If needed, make any adjustments to the positioning.

Explain possible problems to the customer and how to address them.

Conditions of Use

16.3

The **All-Round Lite Joystick** is intended for use as installed by the dealer, in accordance to the installation instructions in this manual.

- The foreseen conditions of use are communicated by the dealer or service technician to the user and/or attendant during the first time use.
- If the usage conditions change significantly, please contact your dealer or a qualified service technician to avoid unintended damage.

Maintenance

17.1

These care and maintenance guidelines will keep the hardware in good condition for a longer period of time and prevent damage:

- Check and tighten all fasteners on a regular basis.
- Repair or replace parts as needed.
- Gently remove dust and dirt from hardware with a soft, damp cloth.



CAUTION

Always check all mounting hardware, making sure each fastener is properly tightened before using the All-Round Lite Joystick.

SAFETY

Replace or repair parts as needed.

NOTICE

Improperly installing the device or altering it in any way will void its warranty.

Cleaning

17.2

- use a soft damp cloth to clean the device.
- Ensure that all cleaners are approved for finished steel, aluminum, plastic, and upholstered surfaces.



CAUTION

Do not immerse in water or use excessive amounts of liquid. Do not apply additional lubricants to the moving parts.



WARNING

As dust and dirt could lead to reduced functionality, all parts of the device should be cleaned on a regular basis.

Monthly Check

17.3

Monthly, or as needed, check to ensure:

- all bolts and screws are still tightened firmly;
- The wiring has sustained no damage; and,
- none of the joystick unit's parts have suffered excessive wear or tear.

A yearly check of the wheelchair and its operating system by a qualified technician is recommended.

| |
|--|
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |

| |
|--|
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |
| |



© 2021, Stealth Products, LLC

Warranty Information



Home Page



Online User Manual



Stealth Products, LLC · info@stealthproducts.com · www.stealthproducts.com
(800) 965-9229 | (512) 715-9995 | 104 John Kelly Drive, Burnet, TX 78611